




NASAL MASK

FITS RIGHT FOR AN EASY NIGHT

The latest nasal mask from ResMed is the only mask with ActiveCell[™] technology. The Activa[™] LT's unique cushion features an "active cell" that expands and contracts to accommodate movement and therapy pressure changes throughout the night. The mask creates a gentle, secure seal and a dynamic defense against mask leak.

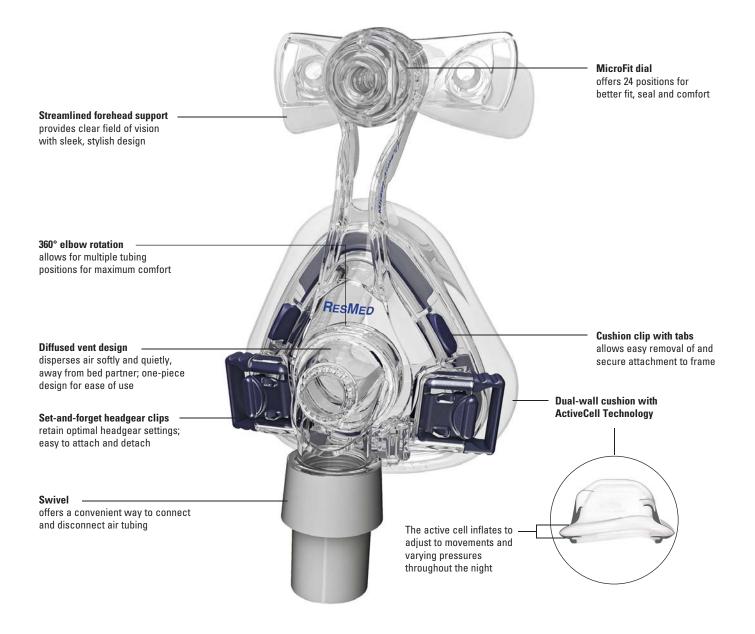

TECHNOLOGY HIGHLIGHT

Activa LT

Same superior Mirage Activa™ sealing technology ... now lighter and more streamlined

Dynamic cushion and MicroFit dial combination ... fitting is quicker and easier than ever

One-piece diffused vent design ... overall quiet performance






A fine-tuned fit

Because the cushion's active cell does all the work, the mask is less dependent on headgear tension to achieve and maintain an effective seal. In fact, the Activa LT works best with some slack in the headgear—no more overtightening straps! That means quicker fitting time, less callbacks and improved user compliance.





With only minor fitting adjustments required, the MicroFit dial provides 24 positions of comfort at your fingertips. Easily adjustable with just one hand (even while lying down), the dial finetunes the fit where users need it most.

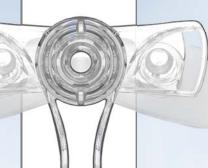





ActiveCell expands ...

comfortable seal.

... and contracts

Room to breathe, freedom to move

ActiveCell for peace of mind The Activa LT cushion is not your average cushion. Once therapy is switched on, the

cushion is drawn to the face, creating a gentle

yet secure seal. As users move or their therapy

pressure changes throughout the night, the cushion's active cell automatically expands and contracts, delivering the most effective,

Because the Activa LT maintains its seal throughout the night, users are free to move around during sleep. The mask automatically, continually and actively adjusts to movement and other changes during sleep, providing a dynamic defense against leak.









Quiet comfort

Featuring an innovative vent design, the Activa LT is not only sophisticated, but also quiet. Located on the elbow for easy cleaning and disassembly, the diffused vents are designed to disperse air quietly and softly, at a short distance, away from the bed partner. The elbow partition also aids CO₂ washout and further reduces noise.





PRODUCT CODES

Mirage Activa LT Nasal Mask

Americas

| Medium | 60148 |
|------------|-------|
| Large | 60149 |
| Large Wide | 60150 |

REIMBURSEMENT CODES

Mirage Activa LT Nasal Mask

| HCPCS | Code Descriptor |
|-------|--|
| A7034 | Nasal mask used with positive airway pressure device, 1 per 3 months |
| A7032 | Replacement cushion for nasal mask, 2 per month |
| A7035 | Headgear used with positive airway pressure device, 1 per 6 month |



ResMed Corp Poway, CA, USA +1 858 746 2400 or 1 800 424 0737 (toll free), ResMed Ltd Bella Vista, NSW, Australia +61 (2) 8884 1000 or 1 800 658 189 (toll free). Offices in Austria, Belgium, Brazil, China, Finland, France, Germany, Greece, Hong Kong, India, Ireland, Italy, Japan, Malaysia, Mexico, Netherlands, New Zealand, Norway, Portugal, Singapore, South Africa, Spain, Sweden, Switzerland, United Kingdom (see website for details). Activa is a trademark of ResMed Ltd. and registered in U.S. Patent and Trademark Office. ©2008 ResMed. 1013193/1 08 09